
EN - English
SV - Svenska
DK - Danske
NO - Norsk
FI - Suomi

Read this Manual carefully before use

SAFETY PRECAUTIONS

- EN
1. People with the following conditions should consult their physician before using the product:
 - Who uses a medical electronic device, which is embedded in the body, such as a pacemaker.
 - Who suffers from heart disease.
 - Who suffers from neck, backbone injury or disease.
 - Who is suffering from bone fragility, especially on the backbone.
 - After surgery.
 - Who is pregnant or shortly after delivery or in a period of monthly illness.
 - Who suffers from malignant tumors.
 - Who is suffering from osteoporosis.
 - Who is suffering from blood circulation problems or thrombosis.
 - Who is suffering from skin disease or infection.
 - Who is undergoing medical treatment, feels some physical abnormality, or is being treated.
 - Who has injuries, diseases due to illness or accident.
 - Who has some abnormality or deformity of limbs or backbone.
 2. Do not let water enter the product or be splashed near the surface of the product.
 3. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or next to swimming pool).
 4. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
 5. If the power cord or product itself is damaged, do not use the product.
 6. To prevent the risk of electric shock, do not open the product outer cover.
 7. Product is not intended for use by the disabled, sensory, or mentally impaired persons (including children), unless there is supervision or guidance by the people responsible for their safety when they use the product.
 8. If the product or the power cord is damaged, do not try to repair the product.
 9. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.
 10. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection.
 11. Do not use the product if cover is damaged or if there is a mechanical or electronic failure.
 12. Only use or place the product in accordance with the instructions in this user manual.
 13. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.
 14. Do not apply medical or cosmetic cream on the product or on your body before using or during the use of this product.
 15. Do not destroy, bend, stretch, twist, or tie the power cord. Damaged power cord can cause fire or electric shock.
 16. Recommended usage time for an individual is 10 minutes at a time. Avoid using the product on the same spot for more than 60 seconds.
 17. Please remove all accessories on body before use.
 18. This product is intended for household use only.
 19. Do not use the product while sleeping, when drunk, or during a fever.
 20. Avoid using the product immediately after eating. Wait for at least an hour before using the product.
 21. Its recommended not to use the product on bare skin to protect muscles if the intensity is too strong.
 22. When massaging, take note of the movement of massage heads and refrain from using excessive intensity. Do not use the product on head, throat, chest, armpit, abdomen, spine, bony areas, joints and inflamed or bruised parts of the body.
 23. Do not use the product while driving or exercising.
 24. If you have long hair, please tie the hair up before using the product.

PRODUCT SUMMARY

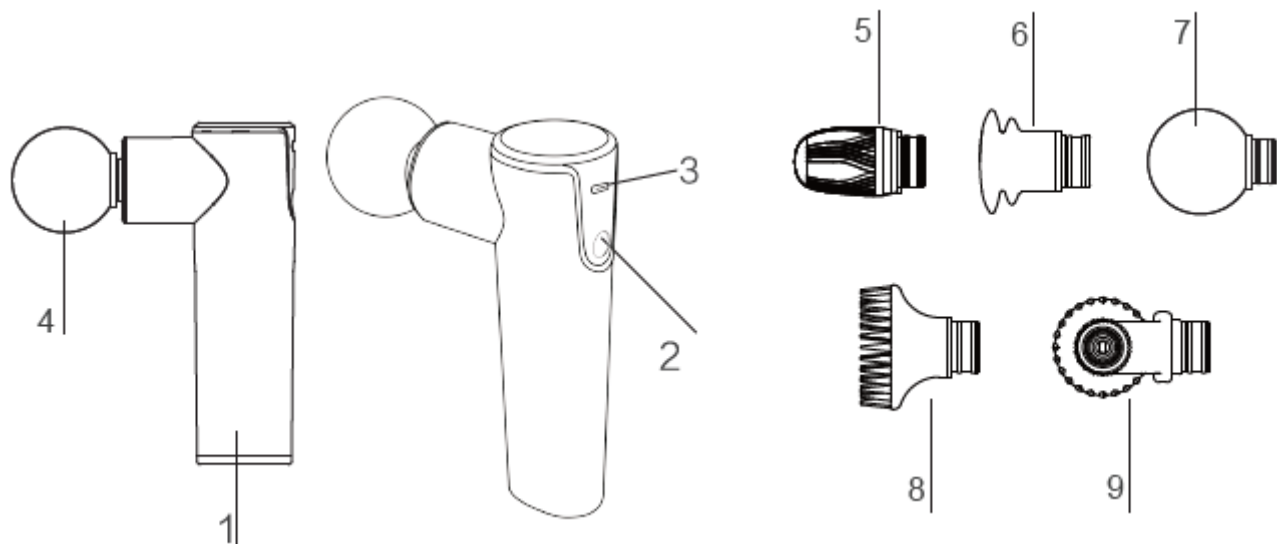1. Main body
2. Power button & Gear selector
3. Charging port
4. Replaceable massage head
5. Bullet Head
6. Pneumatic Head
7. Dome Head
8. Pin-Style Head
9. Rolling Head

DIRECTIONS FOR USE

Caution before use

Before use, please take off your jewelry to avoid discomfort during massage. If you feel uncomfortable during initial use or use, please stop using immediately and consult your physician.

Scope of application

Reduce muscle soreness, promote blood circulation, relieve muscle stiffness, and promote muscle recovery. This product is only suitable for auxiliary relief, if you need treatment, please go to the licensed specialist or hospital for treatment.

About the status after use

The massage area may turn red for a short time, but it will return to normal after a while. The massage area may also leave indentation, but after a period of time will return to normal. When the product is in use, it will stop working or even shut down when it is subjected to excessive resistance. This is the normal motor protection function of the product, and it can be used normally after it is restarted, not because the product is abnormal. It is normal for the main engine to heat up slightly during charging and use.

USE

1. Please charge the product before the first use. While the product is being charged, it cannot be started up for use. Please read chapter "Charging" for more information.
2. For the first time, it is recommended to use it in a low gear to avoid discomfort.
3. Chose massage head, its recommended to use the round or the flat massage for first use. Please see chapter "Change massage head" for how to change.
4. Long press the power button to start the machine. By default, the first gear will be turned on and one indicator will blink.
 - Press to switch the second gear, 1 gear indicator light is steady on;
 - Press again to switch the third gear, and the two gear indicators are always on;
 - Press again to switch the fourth gear, 3 gear indicator light is steady on;
 - Press the fourth time to return to 1, and so on.
5. Press power button for 3 seconds to turn off. The product will automatically shut down after 10 minutes of usage.

Massage levels:

First gear strength: 1 gear indicator light flashing, about 900 revolutions;

Intensity second gear: 1 gear indicator light is always on, the revolution is about 1400;

Intensity third gear: 2 gear indicator lights are always on, the revolution is about 2000;

Intensity fourth gear: 3 gear indicator light is always on, the revolution is about 2800

CHANGING MESSAGE HEAD

1. Remove the massage head by pulling it out. You could feel some vibrations from the massage gun, this is normal. Keep pulling in straight line to remove.
2. Install massage head by inset the grey part of the massage head into the massage guns accessories slot. The grey part of the massage head should not be showed when its correct installed.

FUNCTIONS

Bullet Head - Ideal for muscle release. Designed to target specific smaller muscular knots and focused trigger points.

Pneumatic Head - The soft, squishy feel of the head is designed for managing more painful & sensitive muscles.

Dome Head - General, full body muscle massage & focused muscle work

Pin-Style Head - Soft Massage for Head, Fatigue Relief

Rolling Head - Full body massage by flexible rolling method

CHARGING

1. Plug the original Type-C charging cable into the charging port of the product and connect the other end of the charging cable to the USB port of the 5V-2A adapter.
2. Plug in the charger, the product will start charging, and the three power indicator lights will flash to indicate charging; After charging, the three indicators are always on. After charging, remove the charging device in time.
3. It is recommended to charge the product before the power is used up, around 20-40% is optimal. It's not needed to change to 100%, you can stop at any time. Note: When charging, the machine cannot be turned on. After charging, the charger can be switched on and off normally.
4. To protect bones and muscles and avoid injury caused by long-time fixed-point massage, the intelligent protection mode will be automatically stopped after 10 minutes of continuous use. When gear is switched, the timer will reset.

Power Display:

When 1 light is on, the power is less than 40% (please charge in time);

When 2 lights are on, the power is 40%-70%;

When 3 lights are on, the power is 70%-100%.

STORAGE AND MAINTENANCE

Storage: Store the massage gun in the included storage bag. If stored over 3 months, the machine should be recharged after unpacking and inspection.

Maintenance: Take measures to prevent moisture, rodents, and pollution. The cleaning should be done in the shutdown state with a dry cloth.

PRODUCT PARAMETERS

Name	Massage Gun Compact Power
Model	CHMG400
Output Power	10 W
Rated Time	10 Min
Input Voltage	DC 5V
Battery	1800mA Rechargeable Lithium Battery
Weight	348g
Size	147*75*40mm

Warranty

Champion Nordic guarantees that this product is free from manufacturing defects in terms of materials and workmanship during a period of 2 years in Sweden, Norway, Denmark and Finland. This warranty period begins on the day the product is purchased or delivered. Champion Nordic has no obligation to repair or replace products that are not accompanied by a valid proof of purchase. This warranty applies only to products purchased and used for home use and does not cover damage resulting from abuse, failure to comply with the instructions from Champion Nordic, cases in which the product has been modified or subjected to unauthorized repair, improper packaging on the part of the owner, normal wear and tear, or improper handling on the part of a transport company.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

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Champion ©